**Yes, you can.®**

## *Invacare® Kirilos / Atlante*


---

**User's Manual** **EN**

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**Manuale d'uso** **IT**

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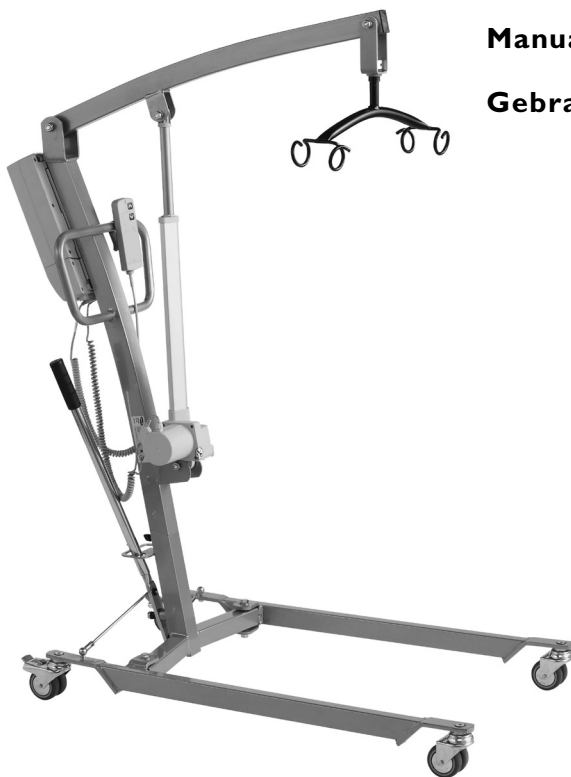*This manual must be given to the user of the product.  
Before using this product, read this manual and save for future reference*

**User's Manual**

**EN**

Stamp of the Distributor

# **Table of Contents**

<b>1.</b>	<b>General</b>	<b>4</b>
	<b>1.1. Symbols</b>	<b>5</b>
<b>2.</b>	<b>Safety</b>	<b>6</b>
<b>3</b>	<b>Technical data</b>	<b>7</b>
<b>4.</b>	<b>Components and functions</b>	<b>10</b>
<b>5.</b>	<b>Set up and installation</b>	<b>11</b>
<b>6.</b>	<b>Using the hoist</b>	<b>12</b>
	<b>6.1. Emergency, lowering and stop</b>	<b>15</b>
<b>7.</b>	<b>Maintenance and cleaning</b>	<b>16</b>
<b>8.</b>	<b>Troubleshooting</b>	<b>18</b>
<b>9.</b>	<b>Product disposal and recycling</b>	<b>19</b>

## Invacare®'s **Kirilos / Atlante**

Thank you for choosing *Invacare®'s Kirilos / Atlante* mobile hoist. The **Kirilos / Atlante** is especially designed to lift the patient to and from wheelchairs, beds and toilets.

*Invacare®* is certified according to DS/EN ISO 9001 and ISO 13485 which ensures that our customers are always supplied with products of uniform quality.

Throughout the entire production process, our materials and components are quality controlled by the operators. If the product does not correspond to the quality demands of *Invacare®*, it will not pass quality control.

In the event of a problem in connection with the delivered product, please contact your local *Invacare®* supplier.

*Invacare®* will take no responsibility if the product is used or assembled in any way other than stated in this user manual. Only accessories mentioned in this user manual may be used on the **Kirilos / Atlante** mobile hoists. This product must only be operated by qualified staff who have received the necessary instruction and training.

Please read this user manual carefully before using the mobile hoist.

## I. General

- **Kirilos / Atlante** is CE-marked in accordance with directive 93/42/EEC concerning medical devices.
- Control unit and motors have been approved according to EN 60601.
- **Kirilos / Atlante** has undergone a risk analysis according to EN/ISO 14971.
- Control unit, motors and hand control are IPX4 protected.
- Lifting capacity: 140 kg (22 stone).
- The hoist has an expected lifetime of 8 years.
- The hoist will not operate when in charge mode.
- Warranty - Frame 3 years, Actuator/ Motor 3 years, Battery 1 year.

If the functions of the hoist change, please see the section on "Maintenance and cleaning". This product must be serviced and maintained by qualified personnel.

If the hoist is used in the area of a shower or bath, ensure that the hoist is wiped clean of any moisture after use. DO NOT store the hoist in a damp area or in a damp condition. Periodically inspect all components of the hoist for signs of corrosion. Replace all parts that are corroded or damaged.

# 1.1. Symbols

The patient is not separated from the ground and the chassis:

Direct current:

Alternating current:

Double insulated:

140 kg (22 stone) max. load (SWL) = (Patient + sling + spreader bar):

The product should be recycled where possible:

Refer to user's manual:

## 2. Safety

Please be aware of the following possible risks when using the hoist:

The hoist must be carefully supervised if the patient is a child.

Please note that it is more difficult to manoeuvre the hoist on carpeted or uneven/textured surfaces as opposed to wooden or flat surfaces; this is due to the starting forces required to initiate the movement that may create additional friction exerted on the castors.

It is recommended not to move the user with the jib at maximum height. Take great care when moving the hoist on carpeted, wet, slippery, rough or uneven surfaces. Never use on slopes.

**Caution!** Never attempt to lift a patient outside the base area of the hoist. When lifting a patient from the floor, the patient must always be located between the two hoist legs, and never outside this area. To avoid injury, utmost care must be taken when raising and lowering a patient.

*Invacare®* accepts no liability for any use, change or assembly of the **Kirilos Atlante** hoist other than as stated in this user manual.

There is a risk of squeezing and entrapment during operation and transport of the hoist. Never exceed the maximum recommended lifting capacity.

It is important that the cables do not become entrapped or jammed in moving parts and castors during operation and transport.


**Important!** *Invacare®* recommend that operation of any mobile lifting equipment is carried out by personnel with prior professional training in moving and handling techniques.

Always choose the sling design and size according to the patient's weight, size and physical ability. Always ensure, that the sling is applied correctly. Be sure to check the sling attachments each time the sling is removed and replaced, to ensure that it is properly attached before the patient is removed from a stationary object (bed, chair or commode). Never leave a patient in the sling unattended. Please ensure that the hoist is not exposed to water, as this can cause damage.

3. Technical data

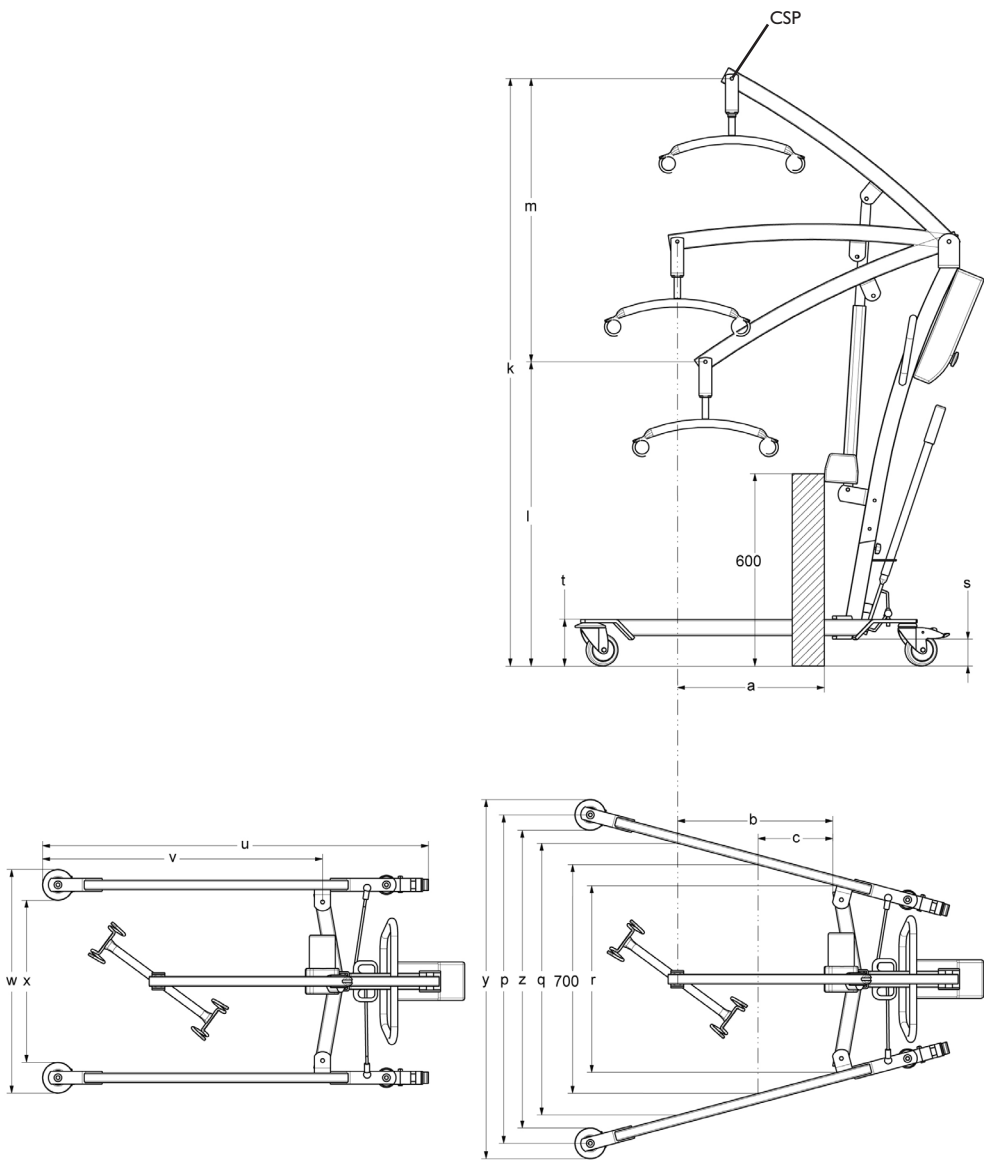General specifications			
	Kirilos	Atlante	
		Electrical	Hydraulic
Highest position (max. height of CSP*) (k)	177 cm	180 cm	
Lowest position (min. height of CSP*) (l)	89 cm	92 cm	
Hoisting range (m)	88 cm		
Lifting area (height range)	60-148 cm	63-151 cm	63-151 cm
Height to upper edge of legs (t)	11 cm	15 cm	
Min. free height (s)	3 cm	6 cm	
Minimum internal width (r)	58 cm		
Internal width at maximum reach (q)	81 cm		
Maximum Length of the base (u)	117 cm	119 cm	
Maximum Internal Length of the Base (v)	86 cm	87 cm	
Maximum External Width of Legs open (y)	103 cm	106 cm	
Total width (open) centre to centre of castors (p)	96 cm		
Minimum Internal Width of Legs (open) (z)	90 cm	87 cm	
Maximum Internal Width of Legs (closed) (w)	70 cm	73 cm	
Minimum Internal Width of Legs (closed) (x)	56 cm	53 cm	
Castor size (FRONT/REAR)	7,5 cm	10 cm	
Sling Material	Polyester	Polyester	Polyester
Maximum lifting capacity (patient+sling+spreader bar)	140 kg (22 stone)		
Total weight (weight out of carton)	35,3 kg	34,1 kg	
Weight, Mast including battery, spreader	19,2 kg	18,1 kg	
Weight, Leg section	16,1 kg		
Maximum reach at 60 cm (a)	48 cm		
Maximum reach from base (b)	50 cm		
Reach from base with legs spread to 70 cm (c)	26 cm		
Turning radius	123 cm		


\* CSP : Central Suspension Point

<b>General specifications</b>	
Operating forces of buttons Max.	Max. 5N
Battery (voltage output)	24V DC === max 8,5A
Charger Input (voltage supply)	230V AC ~ 50 Hz
Charger Output/Charging Time	24V DC === max 800 mA Max 6 hours
Audio/Visual Low Battery Alarm	Yes
Motor Safety Devices	Anti-Entrapment
Working ability	*52 cycles
Velocity of lifting and lowering	<0.15 m/s under maximum load and <0.25 m/s unloaded
Max. current input	Max. 7,0 A DC
Operating temperature	10° C to 40° C
Air humidity	33% to 75%
Sound pressure	< 65 dB
Atmospheric pressure	700 to 1060 hPa
Protection class, control unit	IP20/54
Protection class, hand control	IP66
Protection class, motor	IP66
Insulation class	II - type B
Intermittens	10%, max, 2 minutes/18 minutes
Battery capacity	4,5 Ah
Manual emergency lowering	Yes
Electric emergency lowering/lifting	Yes/No
Expected lifetime	8 years

\*NOTE: Varies depending upon load and stroke.

## 4. Components and functions
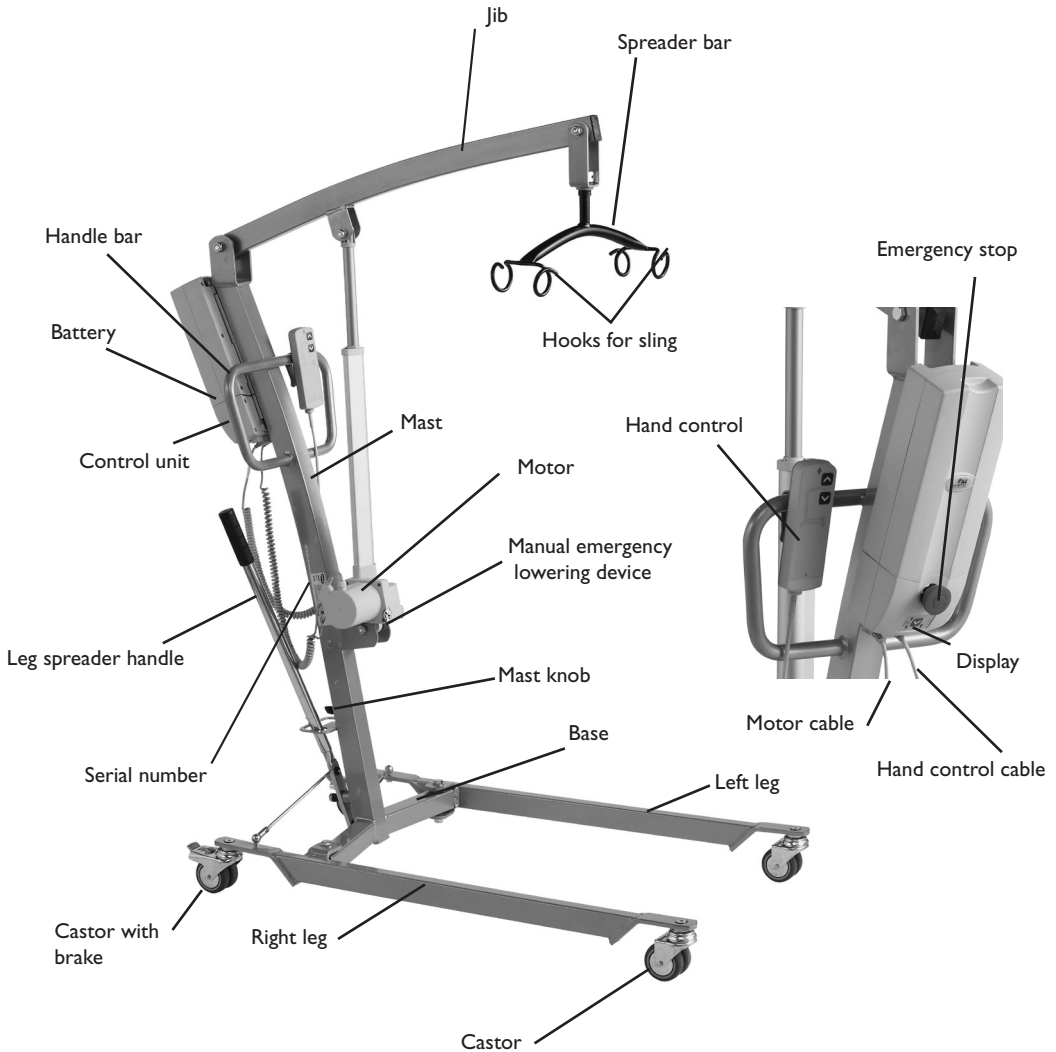
## Accessories

It is compulsory to use original *Invacare*® Accessories which you can purchase through any *Invacare*® dealer.

For repair, please contact your local *Invacare*® dealer.

Service and repair must be carried out by a properly trained technician.

## 5. Set up and installation

### Unpacking

1. Should the packaging be damaged upon receipt, each part of the hoist must be examined for visible defects. In the event of any damage, please contact *Invacare®*.
2. Carefully take out and identify all the parts in the box.
3. The box contains the following parts:
  - 1 hoist
  - 1 mains cable
  - 1 user's manual
  - 1 leg spreader handle
  - 1 battery
  - 1 charger
  - 1 sling

To prevent self-discharging of the battery or accidental operation of the hoist, the emergency stop is activated during shipping. Please remember to deactivate the emergency stop and charge the battery prior to first use.


### Assembling the hoist

The hoist in the Invacare® range are designed for easy assembly. The only tools required are a 17 mm spanner and a 10 mm Allen key

Please adjust the screw of the actuator fixation with the force of 25 Nm.

1. Put the base on a level surface so that all castors make contact with the floor.
2. Lock the rear castors.
3. Loosen the mast knob, but leave it screwed into the base.
4. Insert the tube of the mast into the hole of the base.
5. While supporting the mast assembly, tighten the mast knob.
6. Place the leg spreader handle on its support, push the spring pin to locked the handle on its support.
7. Remove the bolt, assemble the motor and the jib by using the bolt (17 mm key).
8. Plug the hand control cable on the right socket and motor cable on the the left socket of the control unit.
9. Place the battery on the control unit. Push the top of the battery against the mounting bracket until there is an audible click.
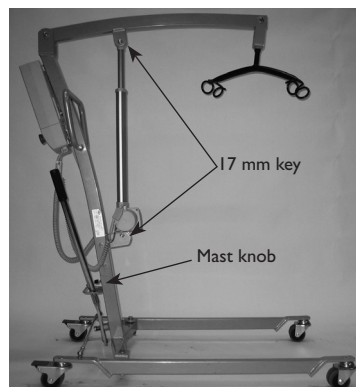
### Dismantling the hoist

1. Put the base on a level surface so that all castors make contact with the floor.
2. Lock the rear castors.
3. Reverse the assembling procedure from steps 9 to 3.

Activate the emergency stop before assembly or disassembly to prevent entrapment/ squeezing. There is a possible risk of squeezing of both limbs and wires during assembly and disassembly of the hoist. Take the utmost care when lifting components during assembly - some parts are heavy. Always remember to adopt the correct lifting position. Make sure there is an audible click when mounting the battery onto the control unit to confirm proper mounting. Otherwise, injury or damage may occur.


## 6. Using the hoist

**Lifting capacity:** 140 kg (22 stone).  
(patient+sling+spreader bar)

### Using the hand control

1. Press “arrow up” on the hand control to raise the jib.
2. Press “arrow down” on the hand control to lower the jib.

**NB** - Operating forces of buttons: Max. 5 N.
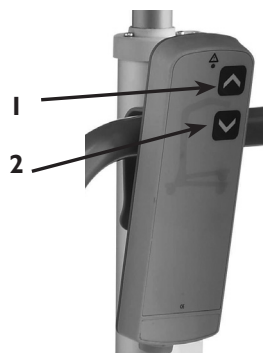
### LED diodes on the control unit

Green :Accumulator full charged

Green/yellow:Accumulator approx.  $\frac{3}{4}$  full

Yellow :Accumulator approx.  $\frac{1}{2}$  full


Yellow/red :Attention! The accumulator should be charged, it is almost empty.

Red :The accumulator is empty and must be charged immediately otherwise a shortening of the service life must be expected. Red light and acoustic signal on pressing the button, it is still possible to perform a complete adjustment cycle. There is no indication when hand control buttons are pressed (except from the remaining battery capacity display or when the emergency lowering/raising system is activated).

Green

Yellow

Red

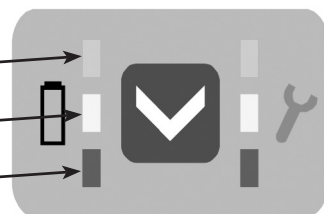
### Operation of the leg spreader handle

The manual leg span is operated by the leg spreader handle on the leg section of the hoist.

The leg spreader handle is operated from one side to the other to adjust the distance between the legs of the hoist.

1. Stand behind the hoist and take a firm hold of the leg spreader handle using both hands.
2. Pull off and turn on the left side to increase the distance between the legs.
3. Turn on the right side to decrease the distance between the legs.
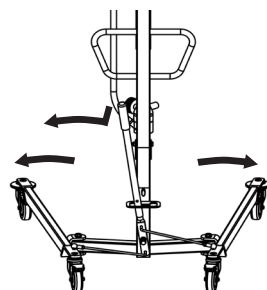
### Brakes

1. Stand behind the hoist and take a good hold of the handle bar.
2. Press the pedals downwards with the foot on the back castors to activate the brake.
3. Press the pedals upwards with the foot to release the brake.

Press down to activate

Push to release

Never try to activate/deactivate the brakes by hand - risk of squeezing !

## Moving the hoist

1. Stand behind the hoist and take a good hold of the handle bar.
2. Release the brakes of the back castors.
3. Now the hoist can be pushed or pulled to the desired location.

## Turning in place


The hoist can be turned (rotated) in place for transfers with limited floor space.

**Note:** To negotiate raised obstacles during movement of the unloaded hoist, it is recommended to pull the hoist backwards.

This operation should only be performed when the hoist is in the unloaded condition.

## Charging the battery

The hoist is equipped with an external charger. It is recommended to charge the batteries regularly to ensure optimal use of the hoist and prolong the life of the batteries. Furthermore, it is recommended to charge the battery before first use.

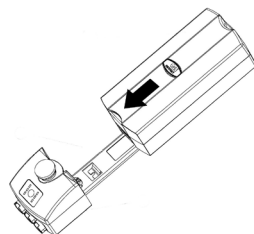The control unit is equipped with a sound signal, which will beep when operating with low battery capacity.

It is recommended to charge the battery as soon as the sound signal is heard.

In addition, **Kirilos / Atlante** is equipped with a display, indicating the approximate remaining battery capacity, when pressing any hand control button.

When installing the battery onto the control unit, it takes a few seconds before the hoist is operational.

## How to charge the battery:

Connect the mains cable to the battery charger and plug it in the socket outlet. It takes about 4 hours to charge the battery. Charging must take place in a room with good air ventilation. The charger stops automatically when the battery are fully charged.

Do not attempt to use the hoist if the battery housing is damaged.  
Replace a damaged battery housing before further use.


## Charging the Battery with battery charger

1. Lift UP on the handle on the back of the battery.
2. Lift the battery up and away from the control box.
3. In order to charge the battery, it be isolated from the control unit. By doing this, the ventilation holes necessary for the charging process are automatically laid bare. Connect the charger and the battery. Finally plug the charger into the house voltage supply. The charge state of battery is indicated by a colored LED.

**NOTE:** The LED lights up constantly green, the battery is charged.

The LED flashes green, the battery is in the process of being charged.

**NOTE:** It will take approximately four hours to charge a battery that requires a full charge.

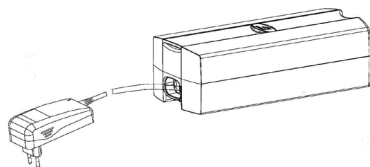4. Disconnect the battery from the battery charger

5. Place the battery on the control unit. Push the top of the battery against the mounting bracket until there is an audible click

**CAUTION: Make sure there is an audible click when mounting battery on the Control Unit. Otherwise, injury or damage may occur.**

### Fitting the sling

Slings designed for 2- or 4-point spreader bars can be used with this mobile hoist, please refer to Invacare®'s sling brochure for details.

**Note:** The size of the sling and the width of the spreader bar should correspond.

The Invacare® Universal Sling can be used for highly physically dependent patients; however, they must have a good degree of head control. The sling allows for easy access for personal hygiene and it is easily applied in both sitting and lying position.

1. Select the sling best suited for the patient in question. Do not attempt to lift any user without prior professional training. If in doubt ask your Invacare® dealer for advice.

2. Fit the sling to the patient as instructed in the Invacare® sling user's guide.

Take great care to ensure that the sling is fitted correctly, so that the patient is comfortable and safe.

3. Manoeuvre the hoist, so that the spreader bar can be lowered to the sling loops as close as possible.

4. Widen the chassis legs for increased stability.

5. Lower the spreader bar; attach the 4 loops of the sling by sliding the loops into the hooks.

6. If the patient is in a sitting position, lower the spreader bar to patient's chest, then move the hoist nearer to the patient. Ensure that brakes are not engaged prior to lifting or lowering.

7. Sitting in chair or wheelchair

Lean the client forward ensuring that they are well supported to reduce risk of falling and to promote feelings of security.

Slide sling down behind the back of the patient in a central position with the aperture apex aligned with the base of the spine. Sit the patient back on against the sling.

Pass the leg support towards the front of the chair to each side of the thigh and under the thigh.

8. Lying in bed, on floor or in bath

Lean the patient forward ensuring that they are well supported to reduce risk of falling and to promote feelings of security.

Slide sling down behind the back of the patient in a central position with the aperture apex aligned with the base of the spine. Sit the patient back on against the sling.


Do not lift or lower a patient with the brakes applied. Always let the hoist find the correct centre of gravity. Never attempt to lift a patient outside the base area of the hoist. When lifting a patient from the floor, the patient must always be located between the two hoist legs, and never outside this area. To avoid injury, utmost care must be taken when raising and lowering a patient. Never exceed the recommended maximum hoist weight and never use the hoist in wet areas.

## Use the correct straps for transfers

**Nota:** make sure that the correct straps are used for transfer


- 1 Black straps: Transfer a patient from the ground
- 2 Green straps: Transfer a patient out of a bed
- 3 Red straps: Other purposes for repositioning a patient
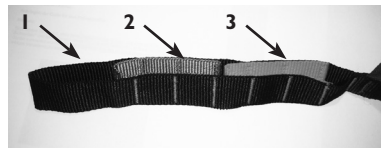
## 6.1. Emergency, lowering, lifting and stop

### Electrical emergency lowering

If the hand control fails, the jib can be lowered by using the circular recessed button for emergency lowering. This is located at the front of the control unit. The jib will be lowered as long as the button is pressed.
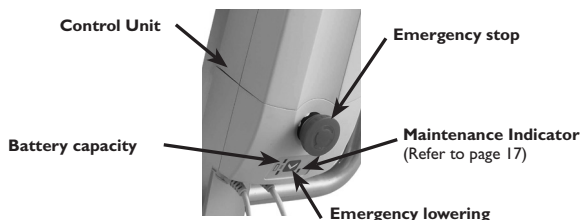
### Manual emergency lowering

In case of partial or total power failure, or if the battery runs down while using the hoist, **Kirilos / Atlante** is equipped with a manual emergency lowering system located at the bottom of the motor.

Turn the 10 mm Allen Key anti clockwise to lowering.

Prior to the manual lowering the key of the power supply of the actuator must be removed.

As long as the actuator is connected with the control box and the power supply, the electrical motor operates as a generative brake.

By interrupting the connection, load and time needed are reduced significantly.


**Note:** Manual emergency lowering is only possible when a patient is sitting in the hoist.

### Emergency stop

In case of emergency: Press the red button for emergency stop.

The emergency stop can be reset by turning the stop button.

10 mm  
Allen key

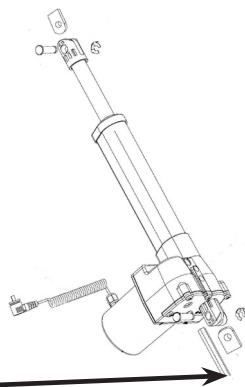**Note!** The mechanical emergency lowering system will only operate when there is a patient in the hoist.

## 7. Maintenance and cleaning

At normal daily operation, a service check-up should take place every year, according to the maintenance chart.

When performing annual or regular maintenance, all parts designed to carry load must be as a minimum tested with maximum load. All safety features must be checked according to EN ISO 10535: 2006 Annex B.

The hoist must be stored at normal room temperature. If it is stored in a damp, cold or wet environment then the motor and other mounting parts may be prone to corrosion.

- The hoist has been designed so that it requires a minimum amount of maintenance.
- It is recommended that the batteries are charged frequently, if possible every night, as it will extend the lifetime of the batteries and will ensure high performance.
- It is recommended that the hoist is cleaned after use.
- Please check that the spreader bar connections, control unit and motors are in good working order and not damaged in any way. If there is any damage to the hoist, refrain from using it, and contact your local *Invacare*® dealer before further use.

### **Wear and damage of loading parts**

It is necessary to maintain and check all parts that are exposed to static or dynamic strain, such as the sling, the spreader bar and the turning points, do not have fractures, are fragile, askew or damaged.

Damaged parts must be replaced before further use.

### **Cleaning**

- To prevent cross-infection, the hoist must be cleaned after use.
- The hoist must be wiped with a moist firmly wrung cloth with ordinary household disinfectants.
- Only use officially approved disinfection detergents.
- Dry the hoist carefully after cleaning.
- Never use acids, alkaline or solvents for cleaning the hoist.

Motors, control unit and mounting parts can be destroyed if the hoist is cleaned in any other way than stated above.

## Maintenance chart

**Only personnel having received the necessary instruction or training by Invacare® must perform service and maintenance on Kirilos / Atlante.**

**Serial number (located on the hoist base):** \_\_\_\_\_

Date:	Initials:					
Visual inspection of all hoist parts						
Inspection of the spreader bar - weldings and assemblies						
Inspection of the motors - undamaged						
Inspection of cables and plugs - undamaged						
Inspection of emergency stop, electrical emergency lowering and raising						
Inspection of the castors - attachment, braking ability and free rolling						

*Invacare®* offers courses in service and maintenance of the hoist.

For all servicing and maintenance requirements, please contact your local *Invacare®* dealer.

## 8. Troubleshooting

Only personnel having received the necessary instruction or training by *Invacare®* must perform service and maintenance on **Kirilos / Atlante**.

Symptom	Possible cause	Remedy
The mast seems loose	The mast screws has become loose	Tighten the screw
Castors are noisy	Dust or dirt in the castors	Clean the castors
Hinge between jib and mast is noisy during operation	Missing lubrication	Lubricate the hinge
The motor is not running	Hand control or motor is not connected	Insert the plugs fully
	No power on the battery	Charge or replace the battery
	The red emergency stop is activated	Turn the red button clockwise to deactivate
Motor noise, but no movement in the piston rod	Motor is damaged	Replace the motor
The piston rod only moves in and not out	Motor is damaged	Replace the motor
The jib motor stops during a patient transfer	Max. load is exceeded	Reduce the load (and the hoist will function normally)

Contact your dealer if the above does not solve your problems.

Lubricate with medically clean oil, e.g. Kemitura Kem Lub KEM-WO 50, order no. 813239.

## 9. Product disposal and recycling

This product has been supplied from an environmentally aware manufacturer that complies with the Waste Electrical and Electronic Equipment (WEEE) Directive 2002/96/CE.

This product may contain substances that could be harmful to the environment if disposed of in places (landfills) that are not appropriate according to legislation.

The 'crossed out wheelie bin' symbol is placed on this product to encourage you to recycle wherever possible.

Please be environmentally responsible and recycle this product through your recycling facility at the end of its life.

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**Invacare® Austria GmbH**
**Invacare® Poirier SAS**
**Invacare® Ireland Ltd**
**Invacare® AB**
**Invacare® B.V.**
**Invacare® Ltd**